


## Anti-Keratin 10 Antibody (KRT10) [clone KRT10/1275] (V3181)

Catalog No.	Formulation	Size
V3181-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	100 ug
V3181-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	20 ug
V3181SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug
V3181IHC-7ML	Prediluted in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide; *For IHC use only*	7 ml

### Bulk quote request

<b>Availability</b>	1-3 business days
<b>Species Reactivity</b>	Human
<b>Format</b>	Purified
<b>Host</b>	Mouse
<b>Clonality</b>	Monoclonal (mouse origin)
<b>Isotype</b>	Mouse IgG1, kappa
<b>Clone Name</b>	KRT10/1275
<b>Purity</b>	Protein G affinity chromatography
<b>UniProt</b>	P13645
<b>Localization</b>	Cytoplasmic
<b>Applications</b>	Immunohistochemistry (FFPE) : 0.1-0.2ug/ml for 30 min at RT
<b>Limitations</b>	This anti-Keratin 10 antibody is available for research use only.

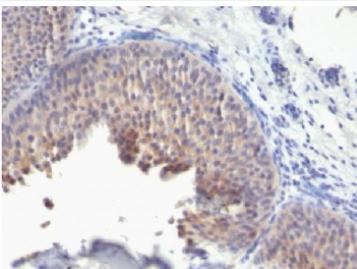

Formalin-fixed, paraffin-embedded human bladder carcinoma stained with anti-Keratin 10 antibody (KRT10/1275).



## Description

This mAb recognizes a protein of 56.5kDa, identified as cytokeratin 10 (CK10). CK10 is expressed in all suprabasal layers of the epidermis. In the epidermis, expression of CK10 strictly parallels the extent of differentiation; it is absent in the basal layer, appears in the first suprabasal layers and increases in concentration towards the granular layer. However, CK10 is rarely detected in early stages of vulvar squamous carcinomas (tumors less than 2 cm, clinical stage I) regardless of the tumor grade. In larger and more advanced tumors (greater than 2 cm, clinical stages II and III), CK10 is detected very frequently. Expression of CK10 is related to maturation of malignant keratinocytes, being preferentially detected in more-differentiated parts.

## Application Notes

Optimal dilution of the anti-Keratin 10 antibody should be determined by the researcher.

1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes
2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

## Immunogen

Recombinant human KRT10 protein was used as the immunogen for the anti-Keratin 10 antibody.

## Storage

Store the anti-Keratin 10 antibody at 2-8°C (with azide) or aliquot and store at -20°C or colder (without azide).