


Anti-Myogenin Antibody [clone SPM144] (V9053)

| Catalog No. | Formulation | Size |
|----------------|---|--------|
| V9053-100UG | 0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide | 100 ug |
| V9053-20UG | 0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide | 20 ug |
| V9053SAF-100UG | 1 mg/ml in 1X PBS; BSA free, sodium azide free | 100 ug |
| V9053IHC-7ML | Prediluted in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide; *For IHC use only* | 7 ml |

Bulk quote request

| Availability | 1-3 business days |
|--------------------|---|
| Species Reactivity | Human |
| Format | Purified |
| Clonality | Monoclonal (mouse origin) |
| Isotype | Mouse IgG1, kappa |
| Clone Name | SPM144 |
| Purity | Protein G affinity chromatography |
| UniProt | P15173 |
| Gene ID | 4656 |
| Localization | Nuclear |
| Applications | Immunohistochemistry (FFPE): 1-2ug/ml for 30 min at RT |
| Limitations | This anti-Myogenin antibody is available for research use only. |

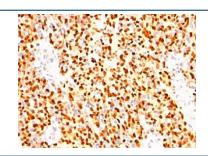

IHC: Formalin-fixed, paraffin-embedded human rhabdomyosarcoma stained with anti-Myogenin antibody (clone SPM144).

Description

Myogenin is a member of the MyoD family of myogenic basic helix-loop-helix (bHLH) transcription factors that also includes MyoD, Myf-5, and MRF4 (also known as herculinor Myf-6). MyoD family members are expressed exclusively in skeletal muscle and play a key role in activating myogenesis by binding to enhancer sequences of muscle-specific genes. The regulatory domain of MyoD is approximately 70 amino acids in length and includes both a basic DNA binding motif and a bHLH dimerization motif. MyoD family members share about 80% amino acid homology in their bHLH motifs. Anti-myogenin labels the nuclei of myoblasts in developing muscle tissue, and is expressed in tumor cell nuclei of rhabdomyosarcoma and some leiomyosarcomas. Positive nuclear staining may occur in Wilms tumor.

Application Notes

The optimal dilution of the anti-Myogenin antibody for each application should be determined by the researcher.

- 1. Staining of formalin-fixed tissues requires boiling tissue sections in pH 9 10mM Tris with 1mM EDTA for 10-20 min followed by cooling at RT for 20 minutes.
- 2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

Immunogen

Human recombinant protein was used as the immunogen for this anti-Myogenin antibody.

Storage

Store the anti-Myogenin antibody at 2-8oC (with azide) or aliquot and store at -20oC or colder (without azide).