


# CEA Antibody / CEACAM5 [clone C66/195] (V2445)


Catalog No.	Formulation	Size
V2445-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	100 ug
V2445-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	20 ug
V2445SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug

# **Bulk quote request**

Availability	1-3 business days
Species Reactivity	Human
Format	Purified
Clonality	Monoclonal (mouse origin)
Isotype	Mouse IgG1, kappa
Clone Name	C66/195
Purity	Protein G affinity chromatography
UniProt	P06731
Localization	Cytoplasmic and luminal surface
Applications	Immunohistochemistry (FFPE): 0.25-0.5ug/ml for 30 min at RT (1) (2)
Limitations	This CEA antibody is available for research use only.

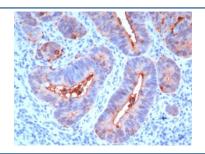

IHC: Formalin-fixed, paraffin-embedded human colon carcinoma stained with CEA antibody (C66/195).

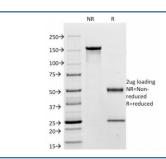

SDS-PAGE Analysis of Purified, BSA-Free CEA Antibody (clone C66/195). Confirmation of Integrity and Purity of the Antibody.

## **Description**

This antibody recognizes proteins of 80-200kDa, identified as different members of CEA family. CEA is synthesized during development in the fetal gut and is re-expressed in increased amounts in intestinal carcinomas and several other tumors. This mAb reacts with nonspecific cross-reacting antigen (NCA) and shows a cross-reaction with human polymorphonuclear leucocytes. It shows no reaction with a variety of normal tissues and is suitable for staining of formalin/paraffin tissues. CEA is not found in benign glands, stroma, or malignant prostatic cells. Antibody to CEA is useful in detecting early foci of gastric carcinoma and in distinguishing pulmonary adenocarcinomas (60-70% are CEA+) from pleural mesotheliomas (rarely or weakly CEA+). Anti-CEA positivity is seen in adenocarcinomas from the lung, colon, stomach, esophagus, pancreas, gallbadder, urachus, salivary gland, ovary, and endocervix.

## **Application Notes**

Optimal dilution of the CEA antibody to be determined by the researcher.

1. Staining of formalin-fixed tissues requires boiling tissue sections in pH 9 10mM Tris with 1mM EDTA for 10-20 min followed by cooling at RT for 20 min

#### **Immunogen**

Recombinant full-length human protein was used as the immunogen for the CEA antibody.

#### **Storage**

Store the CEA antibody at 2-8oC (with azide) or aliquot and store at -20oC or colder (without azide).