


## Chromogranin A Antibody Cocktail [clone LK2H10 + PHE5 + CGA414] (V2098)


Catalog No.	Formulation	Size
V2098-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	100 ug
V2098-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	20 ug
V2098SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug
V2098IHC-7ML	Prediluted in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide; *For IHC use only*	7 ml

### Bulk quote request

Species Reactivity	Human, Mouse, Rat
Format	Purified
Clonality	Monoclonal (mouse origin)
Isotype	Mouse IgG1, kappa
Clone Name	LK2H10 + PHE5 + CGA414
Purity	Protein G affinity chromatography
Buffer	1X PBS, pH 7.4
Gene ID	1113
Localization	Finely granular cytoplasmic
Applications	Immunohistochemistry (FFPE) : 0.1-0.2ug/ml for 30 min at RT
Limitations	This <b>Chromogranin A antibody</b> is available for research use only.

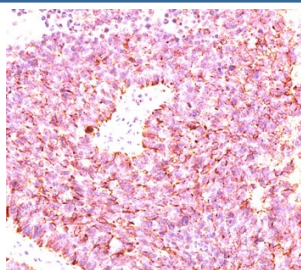

IHC testing of human small cell lung carcinoma stained with chromogranin A antibody cocktail (clones LK2H10 + PHE5 + CGA414).

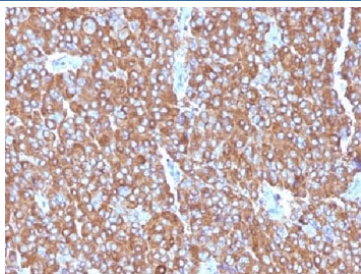

IHC testing of FFPE human adrenal gland with Chromogranin A antibody cocktail.

## Description

Chromogranin A is present in neuroendocrine cells throughout the body, including the neuroendocrine cells of the large and small intestine, adrenal medulla and pancreatic islets. It is an excellent marker for carcinoid tumors, pheochromocytomas, paragangliomas, and other neuroendocrine tumors. Co-expression of chromogranin A and neuron specific enolase (NSE) is common in neuroendocrine neoplasms. Reportedly, co-expression of certain keratins and chromogranin indicates neuroendocrine lineage. The presence of strong anti-chromogranin staining and absence of anti-keratin staining should raise the possibility of paraganglioma. The co-expression of chromogranin and NSE is typical of neuroendocrine neoplasms. Most pituitary adenomas and prolactinomas readily express chromogranin.

## Application Notes

The concentration stated for each application is a general starting point. Variations in protocols, secondaries and substrates may require the antibody to be titrated up or down for optimal performance.

1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes.
2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

## Immunogen

Human pheochromocytoma cells (LK2H10 + PHE5) and recombinant human Chromogranin A protein (CGA414) were used as the immunogen for this antibody.

## Storage

Store the Chromogranin A antibody cocktail at 2-8°C (with azide) or aliquot and store at -20°C or colder (without azide).

## Alternate Names

Beta-Granin; CGA; CHGA; Chromogranin A Parathyroid Secretory Protein 1; ER-37; Pancreastatin; Parastatin; Parathyroid Secretory Protein 1; Pituitary Secretory Protein I; SP-I; Vasostatin I or II, Chromogranin A antibody

## References (2)