


ENO2 Antibody / Neuron Specific Enolase / Gamma Enolase [clone rENO/8684] (V4225)

Catalog No.	Formulation	Size
V4225-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced), 0.05% sodium azide	100 ug
V4225-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced), 0.05% sodium azide	20 ug
V4225SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug

Recombinant **MOUSE MONOCLONAL**

[Bulk quote request](#)

Availability	1-3 business days
Species Reactivity	Human
Format	Purified
Host	Mouse
Clonality	Recombinant Mouse Monoclonal
Isotype	Mouse IgG1, kappa
Clone Name	rENO/8684
Purity	Protein A/G affinity
UniProt	P09104
Localization	Cytoplasm
Applications	Immunohistochemistry (FFPE) : 1-2ug/ml for 30 minutes at RT
Limitations	This ENO2 antibody is available for research use only.

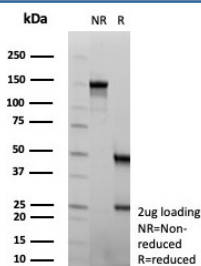

SDS-PAGE analysis of purified, BSA-free ENO2 antibody (clone rENO/8684) as confirmation of integrity and purity.

Description

This monoclonal antibody recognizes a protein of about 50kDa, which is identified as gamma-enolase. Three isoenzymes

of enolases are identified, alpha, beta and gamma. Alpha-isoform is expressed in most tissues, whereas beta-form is expressed predominantly in muscle tissue whereas gamma-enolase is found only in nervous tissue. These isoforms exist as both homodimers and heterodimers, and they play a role in converting phosphoglyceric acid to phosphoenolpyruvic acid in the glycolytic pathway. NSE-gamma is a useful marker to identify peripheral nerves and tumors of neuro-endocrine origins, such as pheochromocytomas. It is usually employed in combination with other markers such as Synaptophysin, Chromogranin A, and Neurofilament.

For detection of neuron-specific enolase (NSE), also known as gamma enolase, across tissue types, see our [NSE antibody](#).

Application Notes

Optimal dilution of the ENO2 antibody should be determined by the researcher.

Immunogen

A recombinant fragment of human NSE gamma (within amino acids 416-433) was used as the immunogen for the ENO2 antibody.

Storage

Aliquot the ENO2 antibody and store frozen at -20oC or colder. Avoid repeated freeze-thaw cycles.