


Kappa Light Chain Antibody [clone KC3L-2] (V3681)

| Catalog No. | Formulation | Size |
|----------------|---|--------|
| V3681-100UG | 0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide | 100 ug |
| V3681-20UG | 0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide | 20 ug |
| V3681SAF-100UG | 1 mg/ml in 1X PBS; BSA free, sodium azide free | 100 ug |
| V3681IHC-7ML | Prediluted in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide; *For IHC use only* | 7 ml |

[Bulk quote request](#)

| | |
|---------------------------|---|
| Species Reactivity | Human |
| Format | Purified |
| Host | Mouse |
| Clonality | Monoclonal (mouse origin) |
| Isotype | Mouse IgG1, kappa |
| Clone Name | KC3L-2 |
| Purity | Protein G affinity chromatography |
| Gene ID | 3514 |
| Localization | Cell surface, cytoplasmic and secreted |
| Applications | Immunohistochemistry (FFPE) : 1-2ug/ml for 30 min at RT Prediluted IHC Only Format : incubate for 30 min at RT (1) |
| Limitations | This Kappa light chain antibody is available for research use only. |

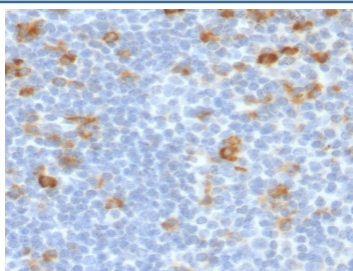

IHC testing of FFPE human tonsil stained with Kappa light chain antibody (clone KC3L-2). Required HIER: boil sections in 10mM citrate buffer, pH6, for 10-20 min followed by cooling at RT for 20 min.

Description

This antibody is specific to the kappa light chain of immunoglobulin and shows no cross-reaction with lambda light chain or any of the five heavy chains. In mammals, the two light chains in an antibody are always identical, with only one type of light chain, kappa or lambda. In general, the ratio of Kappa to Lambda is 3:1. However, with the occurrence of multiple myeloma or other B-cell malignancies this ratio is disturbed. Antibody to the kappa light chain is reportedly useful in the identification of leukemias, plasmacytomas, and certain non-Hodgkin's lymphomas. Demonstration of clonality in lymphoid infiltrates indicates that the infiltrate is malignant.

Application Notes

The concentration stated for each application is a general starting point. Variations in protocols, secondaries and substrates may require the Kappa light chain antibody to be titrated up or down for optimal performance.

1. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

Immunogen

Recombinant human Kappa light chain was used as the immunogen for this Kappa light chain antibody.

Storage

Store the Kappa light chain antibody at 2-8°C (with azide) or aliquot and store at -20°C or colder (without azide).

References (1)