


## Ki-67 Antibody [clone MKI67/8743] (V4109)


| Catalog No.    | Formulation   | Size   |
|----------------|---|--------|
| V4109-100UG    | 0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced), 0.05% sodium azide | 100 ug |
| V4109-20UG     | 0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced), 0.05% sodium azide | 20 ug  |
| V4109SAF-100UG | 1 mg/ml in 1X PBS; BSA free, sodium azide free                          | 100 ug |

[Bulk quote request](#)

|                    |   |
|--------------------|---|
| Availability       | 1-3 business days   |
| Species Reactivity | Human   |
| Format             | Purified  |
| Host               | Mouse   |
| Clonality          | Monoclonal (mouse origin)                                   |
| Isotype            | Mouse IgG   |
| Clone Name         | MKI67/8743  |
| Purity             | Protein A/G affinity  |
| UniProt            | P46013  |
| Localization       | Nucleus   |
| Applications       | Immunohistochemistry (FFPE) : 1-2ug/ml for 30 minutes at RT |
| Limitations        | This Ki-67 antibody is available for research use only.     |

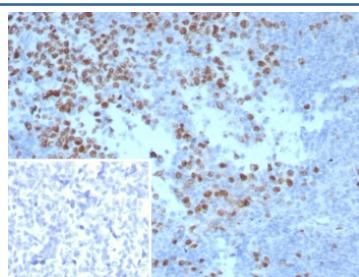

IHC staining of FFPE human tonsil tissue with Ki67 antibody (clone MKI67/8743). Inset: PBS used in place of primary Ab (secondary Ab negative control). HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

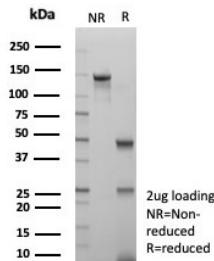

SDS-PAGE analysis of purified, BSA-free Ki67 antibody (clone MKI67/8743) as confirmation of integrity and purity.

## Description

Ki-67 antigen is a nuclear, non-histone protein that is present in all stages of the cell cycle except G0. This characteristic makes Ki-67 an excellent marker for proliferating cells and is commonly used as one of the prognostic factors in cancer studies. A correlation has been demonstrated between Ki-67 index and the histo-pathological grade of neoplasms. Assessment of Ki-67 expression in renal and ureter tumors shows a correlation between tumor proliferation and disease progression, thus making it possible to differentiate high-risk patients. Ki-67 expression may also prove to be important for distinguishing between malignant and benign peripheral nerve sheath tumors. Ki-67 labeling index has been shown to be a prognostic marker in a number of neoplasms including grade II astrocytoma, oligodendrogloma, colon carcinoma, and breast carcinoma. In general, Ki-67 is a good marker of proliferating cell populations.

## Application Notes

Optimal dilution of the Ki-67 antibody should be determined by the researcher.

## Immunogen

A recombinant human Ki67 protein fragment (within amino acids 1-200) was used as the immunogen for the Ki-67 antibody.

## Storage

Aliquot the Ki-67 antibody and store frozen at -20oC or colder. Avoid repeated freeze-thaw cycles.