


Macrophage Antibody [clone LN-5] (V3066)

Catalog No.	Formulation	Size
V3066-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	100 ug
V3066-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	20 ug
V3066SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug

 Citations (13)

[Bulk quote request](#)

Availability	1-3 business days
Species Reactivity	Human
Format	Purified
Host	Mouse
Clonality	Monoclonal (mouse origin)
Isotype	Mouse IgM, kappa
Clone Name	LN-5
Purity	PEG precipitation
UniProt	Not Known
Localization	Cytoplasmic
Applications	Immunohistochemistry (FFPE) : 0.1-0.2ug/ml for 30 min at RT
Limitations	This Macrophage antibody is available for research use only.

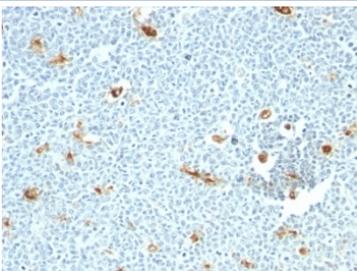

IHC: Formalin-fixed, paraffin-embedded human tonsil stained with Macrophage antibody (LN-5)

Description

This mAb stains the cytoplasm of macrophages and histiocytes in hematopoietic organs, Kupffer s cells of the liver and Langerhan s cells of the skin. Macrophages comprise of many forms of mononuclear phagocytes found in tissues. Mononuclear phagocytes arise from hematopoietic stem cells in the bone marrow. After passing through the monoblast and pro-monocyte states of the monocyte stage, they enter the blood, where they circulate for about 40 hours. They then enter tissues and increase in size, phagocytic activity, and lysosomal enzyme content becoming macrophages. Among the functions of macrophages are nonspecific phagocytosis and pinocytosis, specific phagocytosis of opsonized microorganisms mediated by Fc receptors and complement receptors, killing of ingested microorganisms, digestion and presentation of antigens to T and B lymphocytes, and secretion of a large number of diverse products, including many enzymes including lysozyme and collagenases, several complement components and coagulation factors, some prostaglandins and leukotrienes, and many regulatory molecules (Interferon, Interleukin 1). LN-5 selectively stains human sebaceous glands in formalin-fixed, paraffin-embedded skin samples. Undifferentiated sebocyte progenitors are negative, and only sebocytes from the onset of their differentiation reveal positive cytoplasmic staining. Since there are very few selective and easy-to-use markers of sebaceous glands, LN-5 antibody can offer a simple and relatively specific way to detect human sebocytes from the onset of their.

Application Notes

Optimal dilution of the Macrophage antibody should be determined by the researcher.

1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 min.

Immunogen

Human peripheral blood mononuclear cells were used as the immunogen for the Macrophage antibody.

Storage

Store the Macrophage antibody at 2-8oC (with azide) or aliquot and store at -20oC or colder (without azide).