


TAG-72 Antibody Cocktail [clone B72.3 + CC49] (V3144)


Catalog No.	Formulation	Size
V3144-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	100 ug
V3144-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	20 ug
V3144SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug
V3144IHC-7ML	Prediluted in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide; *For IHC use only*	7 ml

[Bulk quote request](#)

Availability	1-3 business days
Species Reactivity	Human
Format	Purified
Host	Mouse
Clonality	Monoclonal (mouse origin)
Isotype	Mouse IgG1, kappa
Clone Name	B72.3 + CC49
Purity	Protein G affinity chromatography
UniProt	Not Known
Localization	Cytoplasmic and cell surface
Applications	Flow Cytometry : 1-2ug/10 ⁶ cells Immunofluorescence : 1-2ug/ml Immunohistochemistry (FFPE) : 1-2ug/ml for 30 min at RT
Limitations	This TAG-72 antibody cocktail is available for research use only.

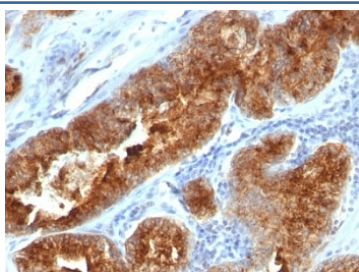

IHC: Formalin-fixed, paraffin-embedded human prostate carcinoma stained with TAG-72 antibody (B72.3 + CC49).

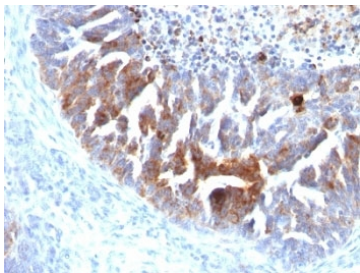

IHC: Formalin-fixed, paraffin-embedded human ovarian carcinoma stained with TAG-72 antibody (B72.3 + CC49).

Description

Recognizes an oncofetal antigen of 220kDa, identified as a tumor-associated glycoprotein (TAG-72) with properties of a mucin. This mAb defines the mucin-carried sialylated-Tn epitope. TAG-72 is usually expressed by adenocarcinomas, but is negative in mesotheliomas. Studies have reported that this antibody has 80% sensitivity and 93% specificity for pulmonary adenocarcinoma. Therefore, TAG-72 is a useful marker to distinguish between mesothelioma and adenocarcinoma. However, false positive reactions can occur so results must be interpreted with the utmost caution. This antibody may be useful in the differentiation of non-small cell carcinomas from small cell carcinomas of the lung. The combined use of anti-TAG-72 and anti-GCDFP-15 is valuable in the diagnosis of apocrine carcinoma.

Application Notes

The optimal dilution of the TAG-72 antibody for each application should be determined by the researcher.

1. Staining of formalin-fixed tissues requires boiling tissue sections in pH 9 10mM Tris with 1mM EDTA for 10-20 min followed by cooling at RT for 20 minutes.
2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

Immunogen

A membrane-enriched fraction of a human breast carcinoma liver metastasis (B72.3) and purified TAG-72 protein (CC49) were used as the immunogens for this TAG-72 antibody cocktail.

Storage

Store the TAG-72 antibody cocktail at 2-8°C (with azide) or aliquot and store at -20°C or colder (without azide).