


# VEGF Antibody / VEGFA [clone VEGF/1063] (V2925)


Catalog No.	Formulation	Size
V2925-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	100 ug
V2925-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	20 ug
V2925SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug

# **Bulk quote request**


Availability	1-3 business days
Species Reactivity	Human
Format	Purified
Clonality	Monoclonal (mouse origin)
Isotype	Mouse IgG1, kappa
Clone Name	VEGF/1063
Purity	Protein G affinity chromatography
UniProt	P15692
Localization	Cytoplasmic, cell surface and extracellular (secreted)
Applications	Immunohistochemistry (FFPE): 1-2ug/ml for 30 min at RT
Limitations	This VEGF antibody is available for research use only.

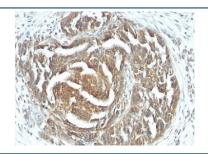


IHC: Formalin-fixed, paraffin-embedded human ovarian carcinoma stained with VEGF antibody (clone VEGF/1063).

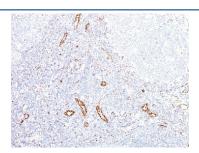

IHC: Formalin-fixed, paraffin-embedded human tonsil stained with VEGF antibody (clone VEGF/1063).

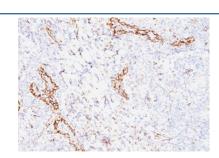

IHC: Formalin-fixed, paraffin-embedded human tonsil stained with VEGF antibody (clone VEGF/1063).

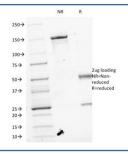

SDS-PAGE Analysis of Purified, BSA-Free VEGF Antibody (clone VEGF/1063). Confirmation of Integrity and Purity of the Antibody.

## **Description**

This mAb recognizes proteins of 19-22kDa (reducing) and 38kDa-44kDa (non-reducing), identified as various isoforms of Vascular Endothelial Growth Factor or Vascular Permeability Factor (VEGF/VPF). It is highly specific to VEGF, which is a homodimeric, disulfide-linked glycoprotein with a close homology to platelet derived growth factor (PDGF). There are multiple isoforms of VEGF containing 206-, 189-, 165-, and 121-amino acid residues. The smaller two isoforms, VEGF165 and VEGF121, are secreted proteins and act as diffusible agents, whereas the larger two remain cell associated. VEGF/VPF plays an important role in angiogenesis, which promotes tumor progression and metastasis.

### **Application Notes**

Optimal dilution of the VEGF antibody should be determined by the researcher.

1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Tris with 1mM EDTA, pH 8, for 10-20 min followed by cooling at RT for 20 min.

#### Immunogen

Human recombinant protein was used as the immunogen for the VEGF antibody.

#### **Storage**

Store the VEGF antibody at 2-8oC (with azide) or aliquot and store at -20oC or colder (without azide).